


## The Next Generation Of Reusable Trocars

Here at Elemis we focus on the next innovations to enhance our trocar devices. These design changes are inspired by physician feedback to create a trocar that will reduce patient trauma and improve the overall procedure.


- Diamond point leading edge to reduce potential patient trauma and soreness
- Diamond point is always sharp (can be re-sharpened) and no longer requires rotation during insertion compared to bevel cut device
- Shorter length to avoid over insertion of pellets past numbed area (prevent the ouch)
- Machined from one solid piece of metal to eliminate bending (others are not reusable past a few times)
- Three piece trocar set that is interchangeable across full line of reusable devices (solving for mismatched sets)
- Increased flat surface area to easily place pellets and limit risk of dropping



**Diamond Point**

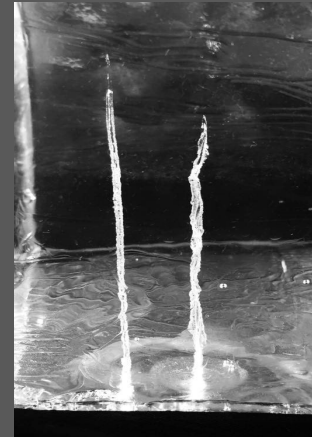

**Bevel Cut**

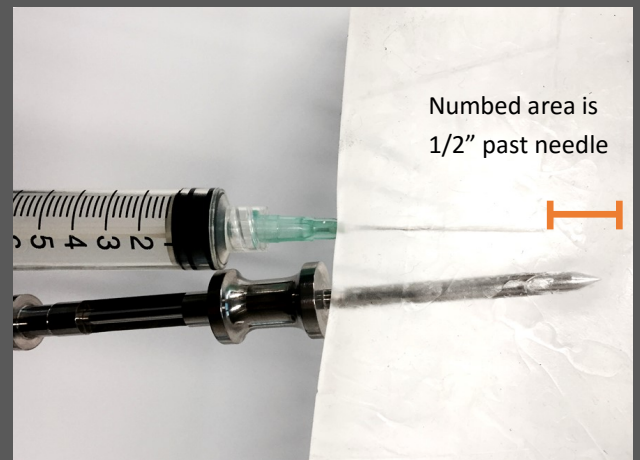

### Elemis Corporation

4010 Technology Way, STE A  
Carson City, NV 89706  
Phone: 775-297-4778  
Fax: 775-297-4776  
Email: [sales@elemiscorp.com](mailto:sales@elemiscorp.com)  
Site: [www.elemiscorp.com](http://www.elemiscorp.com)



Version: 300-03